

Long- & Mixed-Column Nanobore Chromatography for Complex Proteomic Analysis

Christopher J. Toher, Adam W. Perala, Carla Marshall-Waggett, Gary A. Valaskovic

New Objective, Inc., Woburn, MA

Originally published April 2007 at the Association for Biomolecular Resource Facilities Conference, Tampa, Florida, U.S.A.

Introduction

Coupling columns of same or different resin materials is often employed in complex proteomic digest analysis. Despite enhanced separation, these multidimensional columns are costly, time-consuming to produce, and initiate post-column loss by dead-volume introduction. Confounding factors of column-coupling can be eliminated via transparent, true zero-dead-volume (ZDV) unions that achieve flush connections and rapid swap-out facility during system maintenance. In the current investigation, two conventional 10 cm-bedded columns were coupled and connected to the bed terminus of a third 10 cm nanobore column with integrally fritted tip. Analytical merit of this extended column was compared with a single 30 cm-bedded column with integrally fritted tip and the same resin material. These novel unions supported chromatographic data collection with zero dead-volume, negligible resolution loss, and comparable caliber as the single 30 cm-bedded column.

Methods & Materials

Instrumentation & Components

- Ion-trap mass spectrometer (LCQ Deca™, Thermo Fisher Scientific)
- Customized nanospray source (Digital PicoView® 150, New Objective, Inc.)
- NanoLC Pump (Eksigent™)
- Six-port automatic nano-valve (Scivex) with 0.5uL sample loop
- PicoFrit® columns (360 µm OD, 75 µm ID, 15 µm tip ID, New Objective), each containing ProteoPep™ II (New Objective) 5.0 µm-diameter particles packed to 10 cm- and 30 cm- bed lengths
- IntegraFrit™ Columns (360 µm OD, 75 µm ID, New Objective), containing ProteoPep™ II (New Objective) 5.0 µm-diameter particles packed to 10 cm- and 20 cm- bed lengths

Sample Preparation

- A commercially available bovine serum albumin (BSA) standard was diluted to 200 fmol/µL in an aqueous solvent of 2% ACN, 0.1% formic acid
- A commercially available mixture of 5 angiotensins was diluted to 0.1 ng/peptide concentration with 2% ACN, 0.1% formic acid aqueous solvent
- Samples were analyzed via online nanobore ESI-MS in positive-ion-mode

Results

All column combinations were employed in analyzing the angiotensin standard. Data collected using the 30 cm ProteoPep II (PP2)-packed PicoFrit column resulted in FWHMs between 8.4 – 10.2 seconds. The 20 cm IntegraFrit column + 10 cm PicoFrit column combination displayed FWHMs between 13.2 – 14.4 seconds. The two 10 cm IntegraFrit column + 10 cm PicoFrit displayed FWHMs between

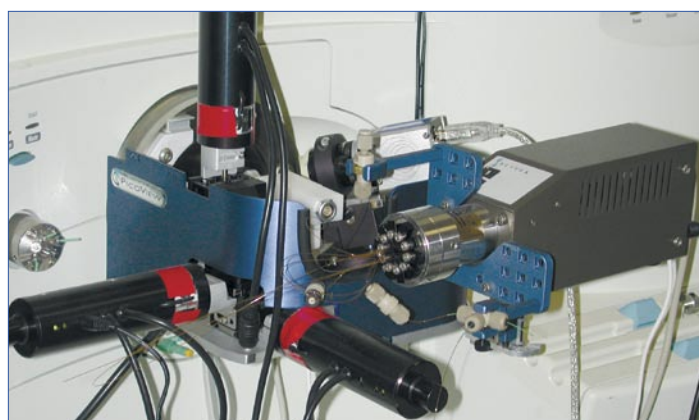


Figure 1 Digital PicoView 150 nanospray source mounted on the Thermo Finnigan LCQ Deca mass spectrometer with Scivex 6-port valve

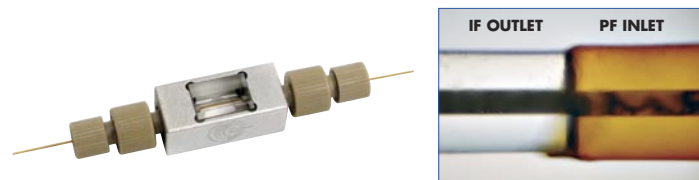


Figure 2 A) PicoClear™ Union, and B) Expanded view of the zero-dead-volume connection achieved inside the clear union body. Note the excellent column-to-column alignment.



Figure 3 Two 10 cm IntegraFrit columns configured with a 10 cm PicoFrit column via two PicoClear Unions to form a single 30 cm column

Angiotensin	MW	Sequence
[Ile ⁷]-Angiotensin III	897.1	RVYIHPI
[Val ⁴]-Angiotensin III	917.1	RVYVHPF
[Asn ¹ ,Val ²]-Angiotensin II	1,031.0	NRVYVHPF
[Val ²]-Angiotensin I	1,282.5	DRVYVHPFHLA
Angiotensin I	1,296.0	DRVYIHPFHL

Table 1 5-Angiotensin composition

12.6 – 14.4 seconds. Figure 5 illustrates three chromatograms from each column combination for analyzing the angiotensin standard; 0.25 ng total peptide were subjected to a 300nL/min flow rate over a 70 minute gradient from 2% - 50% organic modifier concentration.

Figure 6 illustrates the three chromatograms produced in the BSA digest analysis through each column combination; 100fmol BSA were subjected to a gradient identical to that used for angiotensin. Data collected using the 30 cm ProteoPep II (PP2)-packed PicoFrit column allowed 71.8% sequence coverage. The 20 cm IntegraFrit column + 10 cm PicoFrit column supported 58.6% sequence coverage. The two 10 cm IntegraFrit column + 10 cm PicoFrit column yielded 65.1% sequence coverage.

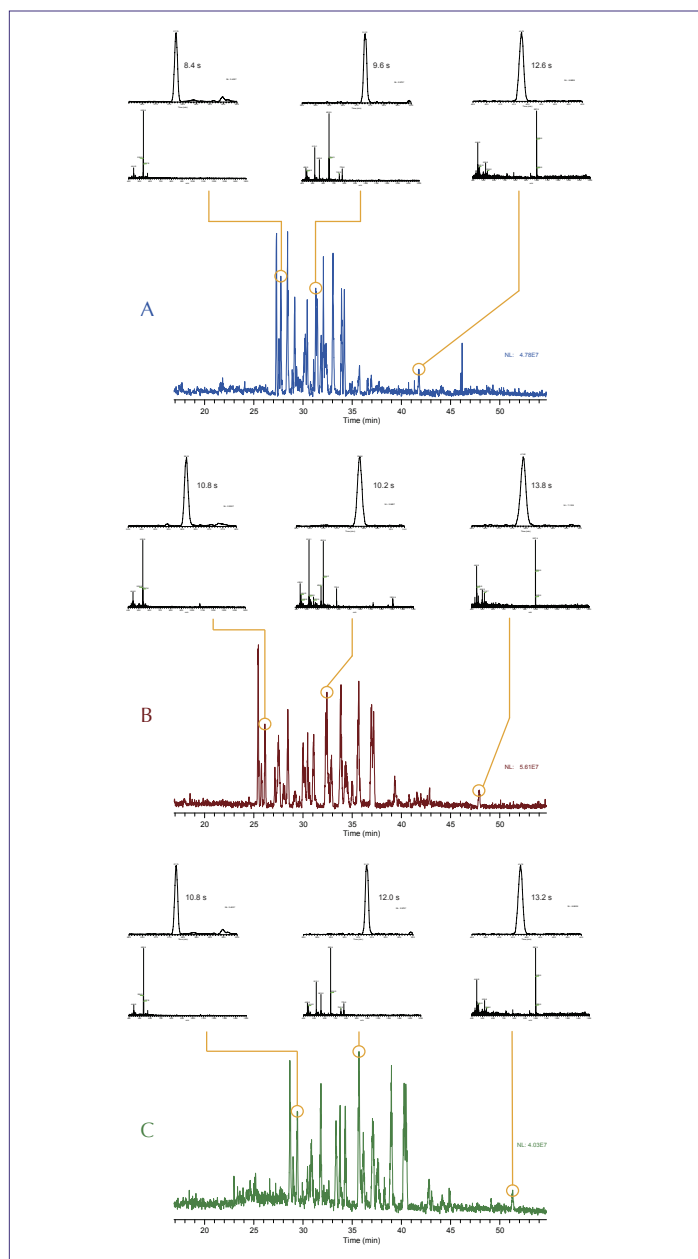


Figure 6 Expanded regions of BSA tryptic digest chromatographic peaks. A) Chromatographic region, as collected with 30 cm PicoFrit column, B) Chromatographic region, as collected with 20 cm IntegraFrit + 10 cm PicoFrit, and C) Chromatographic region, as collected using two 10 cm IntegraFrit columns coupled to a 10 cm PicoFrit column. Injection: 100 fmol BSA, Flow rate: 300 nL/min., Gradient: 2% - 50% B over 70 min.

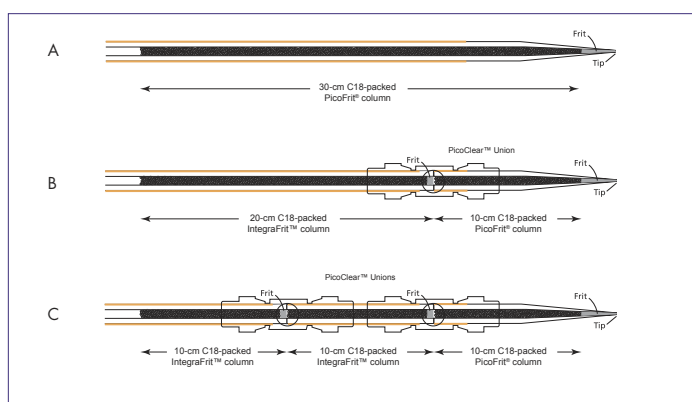


Figure 4 Schematic diagrams of PicoClear union-column combinations. A) PicoFrit column with 30 cm bed, B) 20 cm IntegraFrit coupled to a 10 cm PicoFrit column with a PicoClear Union, and C) Two 10 cm IntegraFrit columns coupled to a 10 cm PicoFrit column via two PicoClear Unions

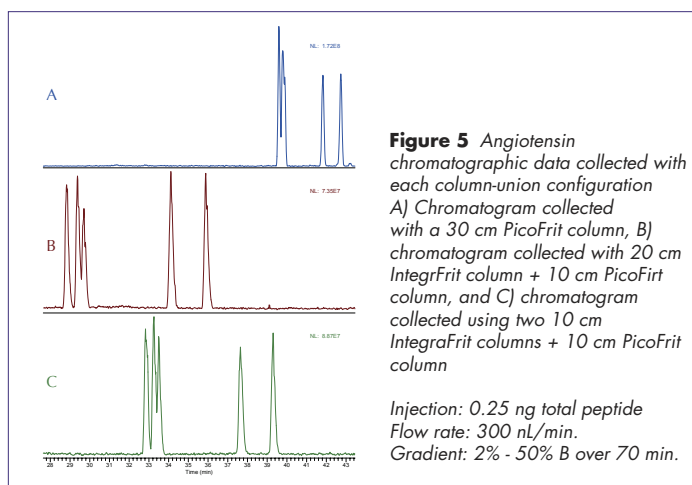


Figure 5 Angiotensin chromatographic data collected with each column-union configuration A) Chromatogram collected with a 30 cm PicoFrit column, B) chromatogram collected with 20 cm IntegraFrit column + 10 cm PicoFrit column, and C) chromatogram collected using two 10 cm IntegraFrit columns + 10 cm PicoFrit column

Injection: 0.25 ng total peptide
Flow rate: 300 nL/min.
Gradient: 2% - 50% B over 70 min.

Conclusions

- Minimal resolution loss and post-column loss were observed for columns combined using transparent, true ZDV unions
- Negligible sequence coverage differences were recorded between each column, although the integral 30 cm column provided the best overall score
- Transparent, true zero-dead-volume (ZDV) unions ensure clean connections between columns without dead volume
- Connecting columns containing different resins will be explored in future work
- Nanobore columns having “semi-disposable” integral guard columns are a viable next step